


How to find the right fit

One guide for all Foal n Pony models — ages 4 to 14+

STEP 1 — UNDERSTAND THE DIMENSION NUMBERS

Every Foal n Pony frame uses a standard three-number code, for example 45 – 16 – 130. Each number tells you one measurement in millimetres. Here is what each one means:

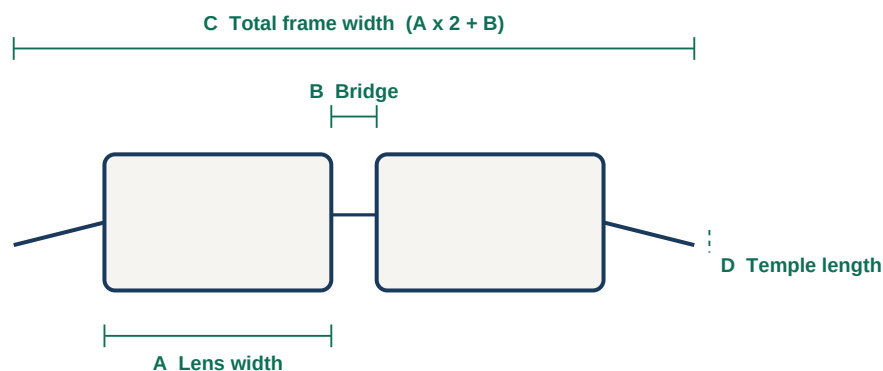

Label	Name	What it means
A	Lens width	Width of one lens opening in mm. E.g. 45mm
B	Bridge	Distance between the two lenses in mm. E.g. 16mm
C	Total width	Lens x 2 + Bridge. Compare to child's face width.
D	Temple	Length of the arm from hinge to tip in mm. E.g. 130mm

STEP 2 — READING A SIZE CODE: EXAMPLE

45 – 16 – 130

Lens: 45mm

Bridge: 16mm

Temple: 130mm

Total frame width = 45 + 45 + 16 = 106mm → measure your child's face and compare to this number.

STEP 3 — FIND YOUR CHILD'S AGE GROUP

Each Foal n Pony model is sized for a specific age group. Use the chart below to find the right range, then check the frame's dimension code.

The most reliable check is always the total frame width (C) vs your child's face width measured temple-to-temple.

YOUNG CHILDREN	Ages 4 – 7			
	YOUNG CHILDREN			
	A Lens width	B Bridge	C Total width	D Temple
	44 – 46 mm	15 – 18 mm	103 – 110 mm	130 – 136 mm
SCHOOL AGE	Head circumference (standard average):			
	49 – 52 cm			
	Models: LUNA · SCOUT · RIVER · JUMPER · PIPPIN · DAISY · FERN			
	<i>Smallest frames in the range. Look for Special Nose Pad or Flexible Side models.</i>			
SCHOOL AGE	Ages 8 – 12			
	SCHOOL AGE			
	A Lens width	B Bridge	C Total width	D Temple
	46 – 48 mm	15 – 18 mm	107 – 114 mm	136 – 137 mm
TWEENS	Head circumference (standard average):			
	52 – 55 cm			
	Models: FABLE · OLIVER · HARPER · PIXIE · CLOVER · BECKETT · SAWYER · SKIPPER · ARCHER · STAR · DAISY · FERN · WILLOW			
	<i>Largest selection. Soft Flex, Medium Flex, and Polarised Clip-on all available.</i>			
TWEENS	Ages 13+			
	TWEENS			
	A Lens width	B Bridge	C Total width	D Temple
	49 – 50 mm	16 – 17 mm	114 – 117 mm	136 – 138 mm
	Head circumference (standard average):			
	54 – 57 cm			
	Models: ARCHER · STAR · ELLE · WILLOW			
	<i>Largest frames. Closer to adult sizing. ARCHER, STAR and WILLOW overlap with 8–12.</i>			

STEP 4 — HOW TO MEASURE YOUR CHILD'S FACE

- 1 Face width**
Use a soft measuring tape. Measure temple to temple across the widest part of the face. This should be close to — or very slightly less than — the total frame width (C).
- 2 Bridge / PD**
Measure the pupillary distance (PD): the space between the centres of both pupils. This should be close to the bridge measurement (B) on the frame.
- 3 Temple fit**
The temple arm (D) should reach comfortably behind the ear without pressing into the skull. If the tip pokes out beyond the ear, the temple is too long.

STEP 5 — CHOOSE THE RIGHT TECHNOLOGY FOR YOUR CHILD

-  **Special Nose Pad**
Ages 4–7 only · LUNA, SCOUT
Best for children wearing glasses for the first time. The special nose pad grips gently without pinching, keeping frames in position even during nap time or active play. Ideal for flat or low nose bridges.
-  **Flexible Side**
Ages 4–14+ · 12 models
Temples bend beyond 90° and return to shape without memory loss. Soft Flex (■400) is slightly more pliable; Medium Flex (■375) has a firmer hold. Best all-round choice for active children.
-  **Polarised Clip-on**
Ages 4–14+ · PIPPIN, DAISY, SAWYER, SKIPPER, WILLOW
Prescription frame with a detachable polarised polycarbonate clip-on. Converts to sunglasses in seconds. Best for children who spend time outdoors or in bright sunlight.

QUICK FIT TIPS

-  The frame should sit level — not tilting up or down on one side.
-  The temples should rest lightly on the ears, not dig in or sit above them.
-  The nose pad (or bridge) should sit flush with no gaps or pressure marks.
-  There should be roughly 5mm clearance between the lash line and the lens.
-  When in doubt, visit a Foal n Pony stockist for a professional fitting.

Still not sure? We're here to help.

WhatsApp us or visit a Foal n Pony stockist near you for a professional fitting.